

Milliken

# PERMEADERM

Translucent Bioactive Matrix


## GLOVE

PRODUCT INFORMATION


### SIZING INSTRUCTIONS

1. Measure from thumb web space (A) to pinkie web space (B) in cm.
2. Add 2 cm.
3. Use chart below to determine correct glove size.

*Example: Thumb web space to pinkie web space = 7.5 cm  
 7.5 cm + 2 cm = 9.5 cm  
 Correct glove size = Medium*

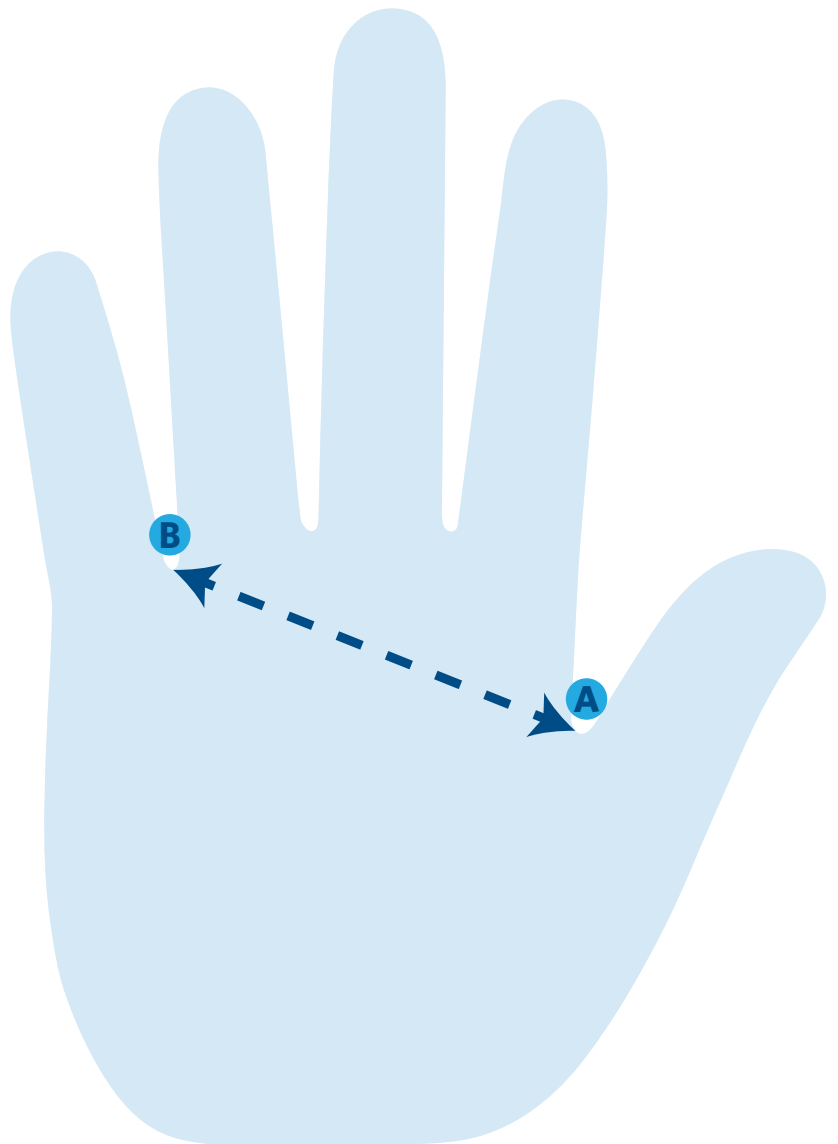

Measurement	+ 2 cm	=	Glove Size
5.5 cm	+2 cm		Extra Small / 7.5 cm
6.5 cm	+2 cm		Small / 8.5 cm
7.5 cm	+2 cm		Medium / 9.5 cm
8.5 cm	+2 cm		Large / 10.5 cm
10.5 cm	+2 cm		Extra Large / 12.5 cm

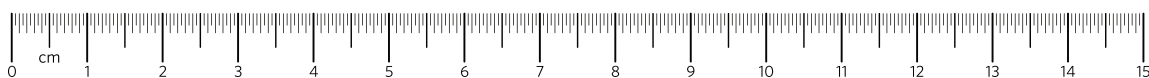

## INDICATIONS FOR USE

### Permeaderm Sheets:

- Partial thickness burn wounds
- Donor sites
- Coverage of meshed autografts

### Permeaderm Gloves:

- Coverage of debrided partial thickness hand burns

## FEATURES & BENEFITS

- Adherence - flexible silicone-coated bioactive matrix allows for optimal wound contact and early mobilization once adherence is achieved
- Transparency - allows for visual wound assessment
- Moisture management - porosity can be varied by applying different amounts of stretch
- One-time application - reduces OR time and the need for multiple dressing changes
- Unique sizes - Multiple glove sizes and large sheets available

## ORDERING INFORMATION

Size	Area	Code	Packaging
PermeaDerm 5" x 10" (12.7 cm x 25.4 cm)	322.6 cm <sup>2</sup>	3000124725	8 ea/box
PermeaDerm 10" x 15" (25.4 cm x 38.1 cm)	967.7 cm <sup>2</sup>	3000124724	3 ea/box
PermeaDerm 15" x 30" (38.1cm x 76.2 cm)	2,903.2 cm <sup>2</sup>	3000124726	1 ea/box
PermeaDerm Glove Extra Small (7.5 cm)	N/A	3000124727	2 ea/box
PermeaDerm Glove Small (8.5 cm)	N/A	3000124728	2 ea/box
PermeaDerm Glove Medium (9.5 cm)	N/A	3000124729	2 ea/box
PermeaDerm Glove Large (10.5 cm)	N/A	3000124730	2 ea/box
PermeaDerm Glove Extra Large (12.5 cm)	N/A	3000124741	2 ea/box

